


Isolator[®] Synergy[™] Access[®]

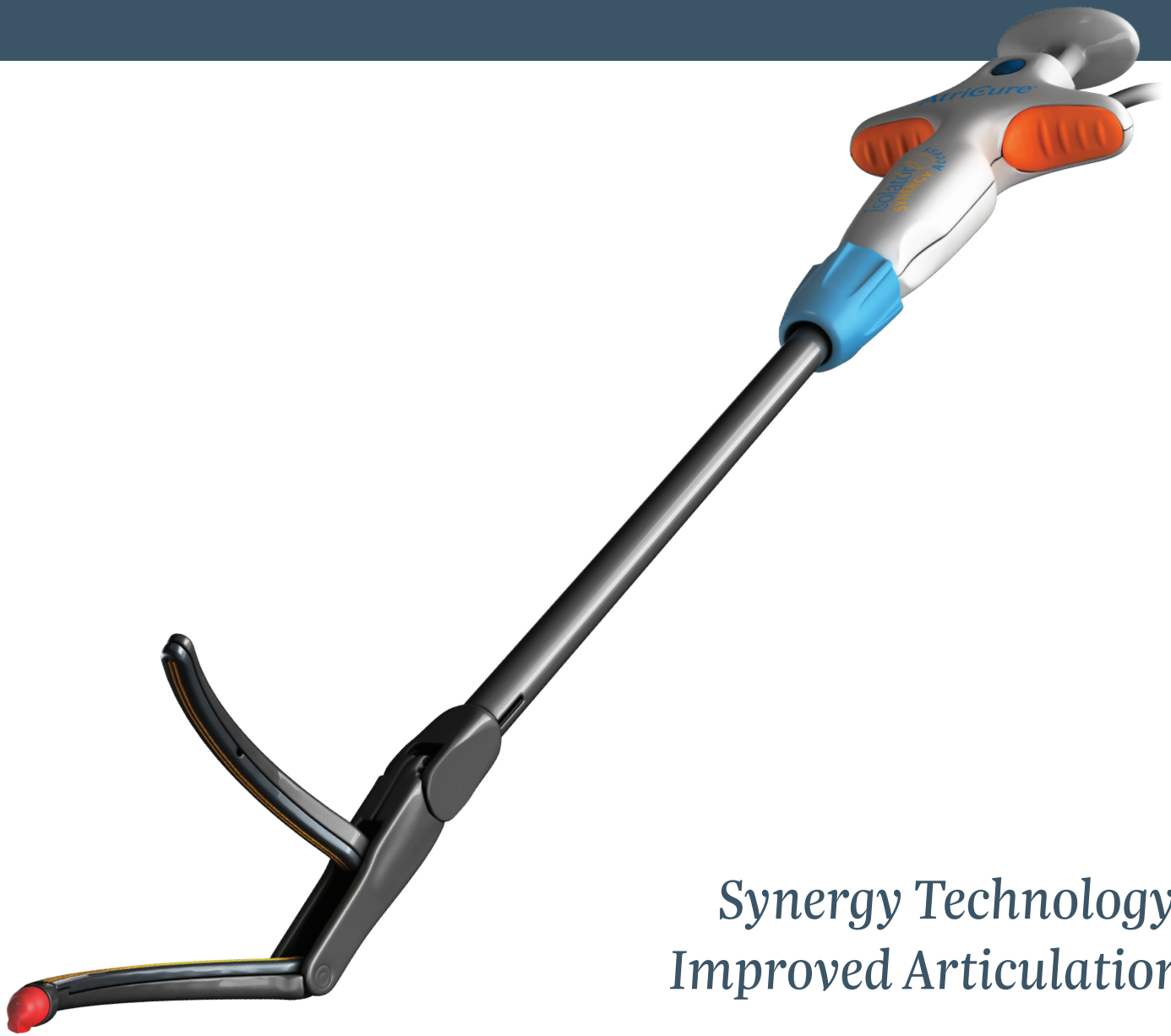


*Synergy Technology,
Improved Articulation*

AtriCure

Isolator Synergy Access

Isolator Synergy Science

Ground-breaking Synergy Science uses two pair of alternating bipolar electrodes to maintain resistive heating between the four poles. An initial lesion is formed where the RF energy fields overlap and continues until a robust, column shaped lesion is produced.

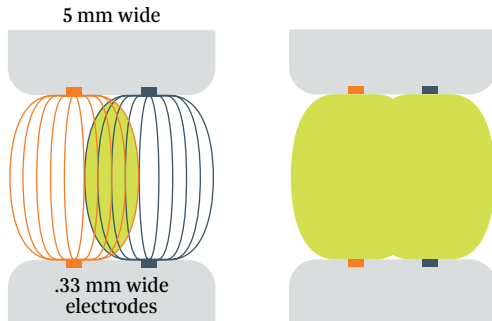

Consistent Tissue Compression

Minimal Tip Splay (<0.06")


+ 0% Jaw Deflection

= True Parallel Closure



Jaw Articulation

The Synergy Access jaw rotates 30 degrees upward and 30 degrees downward to enable ease of access in even the most difficult anatomy.

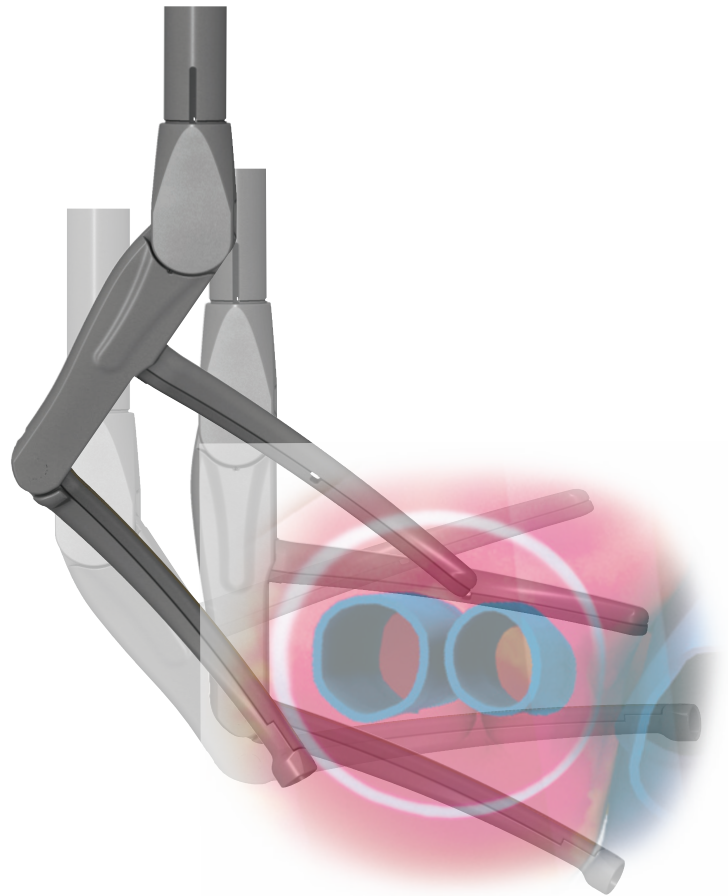

Isolator Synergy Access

Product Code

EMT1

*Includes Glidepath Tape

ATRICURE, INC.

7555 Innovation Way
Mason, OH 45040 USA

+1 (513) 755-4100

+1 (888) 347-6403

CustomerService@AtriCure.com

www.AtriCure.com

PM-US-0114C-0125-G

U.S. Indications: The AtriCure Bipolar (Transpolar) System is intended to ablate cardiac tissue during surgery.

Please review the Instructions for Use for a complete listing of contraindications, warnings, precautions and potential adverse events prior to using these devices.

Rx Only.

AtriCure