


*cryo*SPHERE<sup>®</sup>+




Precision.  
Protection.  
Performance.

AtriCure

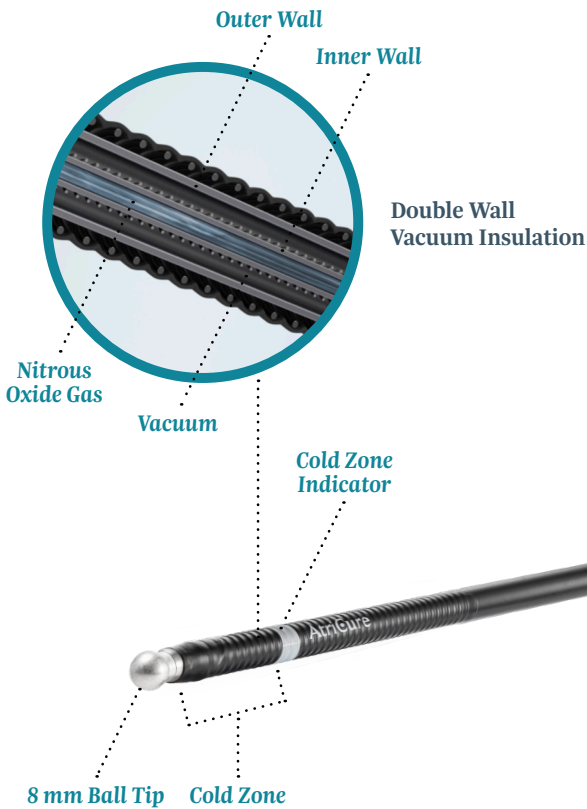

### Cryoablation Probes

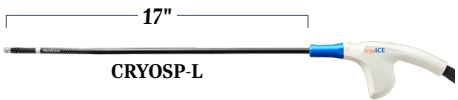

### cryoICE BOX® - V6


ACM (US)

### cryoICE Cart & Upgrade Kits

ASCS (cart with double tank holder)

20 lbs N<sub>2</sub>O tanks not provided by AtriCure

# Cryo Nerve Block: A non-opioid solution for post-operative pain

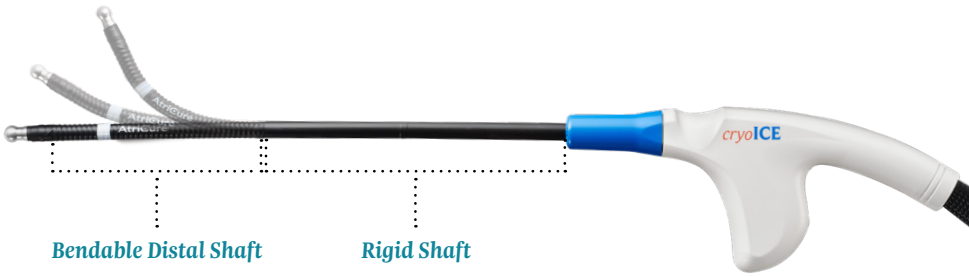

## Features

<b>Double Wall Vacuum Insulation</b>	Designed to reduce risk of injury to non-target tissues that may contact the insulated shaft
<b>Time To Therapeutic Temperature</b>	Optimized time* to reach therapeutic temperature
<b>Cold Zone Indicator</b>	Provides visual reference that is approximately 4 cm* from the distal tip of the probe shaft**
<b>8 mm Ball Tip</b>	8 mm ball tip specifically engineered for cryoablation of peripheral nerves Purposefully shaped to reach target tissue from various approaches/angles
<b>Active Defrost</b>	Allows the probe to be easily removed from the cryoablation site
<b>Rigid Shaft</b>	Stainless steel shaft provides rigidity and tactile feedback for the user
<b>Bendable Distal Shaft</b>	Distal shaft bendable up to 140 degrees

\*Reference on file at AtriCure

\*\*Shaft distal of the Cold Zone Indicator can cause unintended cryoadhesion if touching unintended structures

**U.S. Indications: For Adult Patients** AtriCure's cryoICE® cryoSPHERE® + cryoablation probes are sterile, single use devices intended for use performed by freezing target tissues, creating an inflammatory response (cryonecrosis) for blocking pain by temporarily ablating peripheral nerves. **For Adolescent Patients** The cryoICE cryoSPHERE+ cryo-ablation probes are intended for use to temporarily block pain by ablating intercostal nerves under direct visualization<sup>1</sup> in adolescent patients of at least 12 years of age.

<sup>1</sup>Direct visualization, in this context, requires that the surgeon is able to see the targeted tissue for cryoablation directly or with assistance from a camera, endoscope or other similar optical technology.

**ATRICURE, INC.**

7555 Innovation Way  
Mason, OH 45040 USA

+1 (513) 755-4100

+1 (888) 347-6403

CustomerService@AtriCure.com

www.AtriCure.com

PM-US-3132B-0326-G

**AtriCure**