

*cryo***XT**<sup>™</sup>




Taking on Post-  
Amputation Pain

**AtriCure**

## **cryoXT harnesses cryoICE® technology intraoperatively to manage post-operative pain with multi-surface freezing technology**


cryoXT is designed for use intraoperatively during amputation procedures, featuring a prong-shaped distal tip sized for the nerves found in those applications. Optimal cryoablation is achieved by positioning the nerve entirely within and in contact with the prongs, optimizing heat extraction.

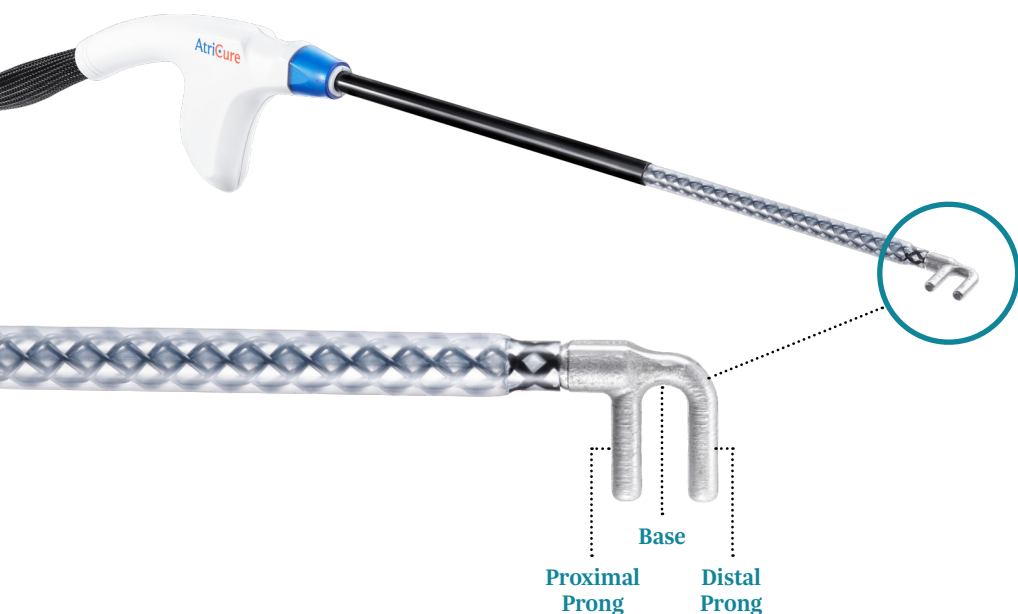

### **Multi-Surface Freezing Technology**

The prong-shaped distal tip is designed to engage exposed peripheral nerves, such as those in amputation procedures.

- Atraumatic fixed, large-bore diameter prongs
- Sized to accommodate nerves found in the extremities
- Surrounds nerve with three freezing surfaces, providing consistent, complete cryoablation

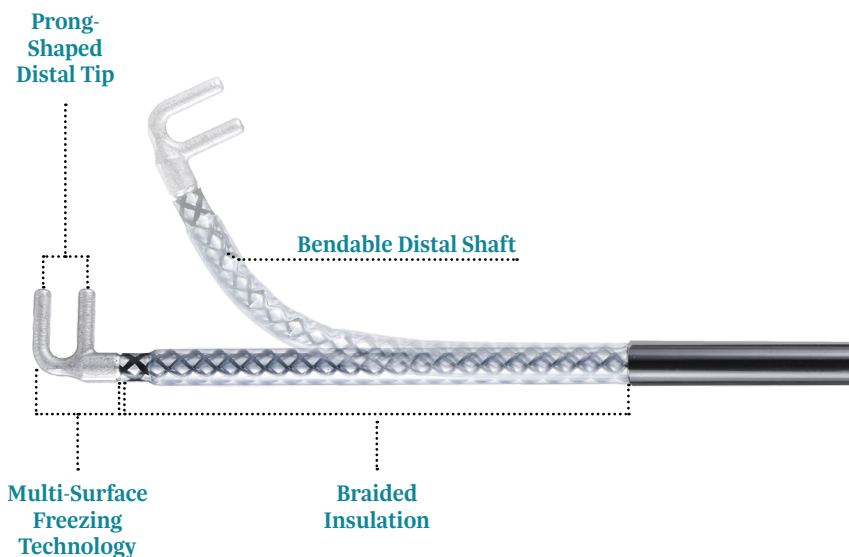

## Features

### Prong-Shaped Distal Tip

- Atraumatic fixed, large-bore diameter prongs
- Designed to engage exposed peripheral nerves
- Sized to accommodate peripheral nerves in extremities

### Multi-Surface Freezing Technology

- Surrounds nerve with three unique surfaces, optimizing heat extraction

### Active Defrost

- Allows probe to be safely and easily removed from cryoablation site

### Bendable Distal Shaft

- Bendable distal shaft up to 90 degrees

### Braided Insulation

- Technology that is designed to maintain a temperature warmer than  $-10^{\circ}\text{C}$  in tissue

### Adaptable

- Multi-functional design allows flexibility in use to accommodate variety of procedures

## Cryoablation Probe



## cryoICE BOX® - V6

ACM (US)

## cryoICE Cart & Upgrade Kits

ASCS (cart with double tank holder)

20 lbs N<sub>2</sub>O tanks not provided by AtriCure

**U.S. Indications:** AtriCure's cryoICE cryoXT cryoablation probes are intended for use to temporarily block pain by ablating peripheral nerves performed by freezing target tissues, creating an inflammatory response (cryonecrosis).

### ATRICURE, INC.

7555 Innovation Way  
Mason, OH 45040 USA

+1 (513) 755-4100

+1 (866) 349-2342

CustomerService@AtriCure.com

www.AtriCure.com

PM-US-4538A-0427-G

**AtriCure**